

CERVMAX™

*A new family of collars...
because each of us is unique.*



Cervmax™ Cervical Orthosis Collar

- Increased comfort promotes patient compliance.
- Provides post-operative and trauma support.
- Engineered to provide cervical support, protection and comfort while limiting range of motion to allow patients to heal properly.
- Multiple sizes to fit a newborn infant to a mature adult.
- Patented color-specific size coding for ease of fitting and more aesthetic appeal to children.

PMT's cervical orthosis collar was engineered with your patients' comfort and stabilization in mind.



Cervmax™ cervical orthosis collars in adult (top) and youth (bottom) sizes.

PMT®

Quality, Innovation, Customer Driven

Product Information



Description:

PMT offers several different models to fit every physician's preference:

- Model 1555** Firm and rigid using a thick plastic shell.
- Model 1565** Soft and comfortable using a thin plastic shell with anterior and posterior reinforcements.
- Model 1575** Firm and rigid using a thick plastic shell with modified shoulder piece for patient comfort.
- Pediatric Model** Specifically engineered for a child's anatomy, not simply a reduced adult model. Patented color-specific size coding allows quick size identification.

Sizes

Adult



Pediatric




The unique V-Back design of the Cervmax™ Collar provides maximum occipital support, cradling the back of the head, and is flexible enough to conform to most neck shapes. Cervmax™ features less shell coverage in the chin area which reduces skin irritation and maintains the maximum mandibular support to control tilt (coronal flexion and extension). The slotted front portion aids in even weight distribution across the sternal notch area, increasing patient comfort with less pressure. Lightweight materials are MRI compatible, silicone-free and latex-free.

How Supplied

Each Cervmax™ Collar is provided standard with an extra set of replacement pads and size identification card.

Our Commitment...Customer Satisfaction

PMT is committed to providing you the highest quality products and the finest customer service. Local sales representatives are available to serve you. Please contact PMT at 1-800-MANKIND (626-5463).

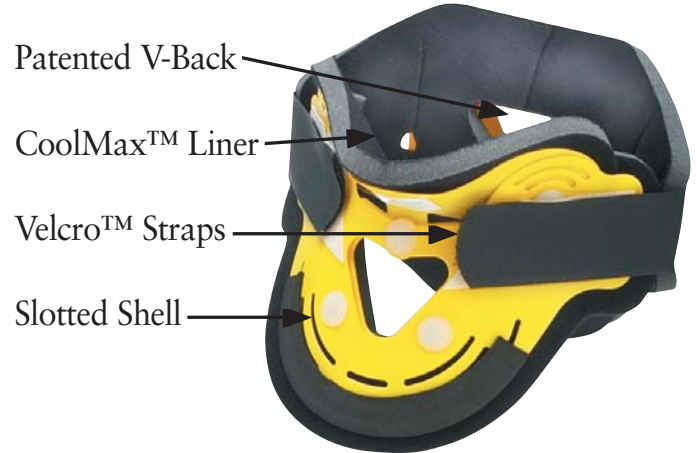

Cervmax™ Features	Cervmax™ Benefits
A. Patented V-Back Design	Provides superior occipital support and immobilization.
B. Patented color-specific size coding	For ease of fitting and more aesthetic appeal to children.
C. CoolMax™ Liner	Wicks moisture away from skin, allows skin to breathe. A hook and loop system allows for easy removal of washable, fast-drying, durable liner.
D. Velcro™ Straps	For easy application.
E. Slotted Shell	Allows airflow and aids in breathability.

PMT and Mankind are registered trademarks and Cervmax™ is a trademark of PMT Corporation. Several Patents Pending. CoolMax™ is a DuPont trademark.
© 2001 PMT Corporation. All rights reserved.