

When you think safety, think PMT Traction Tong



Titanium, sterile
spring loaded
pressure pin

- Strong
- Sterile
- Easy to use

**MRI and CT
compatible**

GraphComp™ 1210 Traction Tong



PMT®
CORPORATION

Quality, Innovation and Customer Driven

1-800-MANKIND



Product Information

Description:

The Model 1210 GraphComp™ Traction Tong puts more diagnostic tools at your disposal including clear MRI and CT scanning. The tong is made of carbon-graphite and the pressure pins are made of titanium, making the set not only MRI and CT compatible, but also exceptionally strong and lightweight for the maximum strength to weight ratio. Model 1210 is available in three sizes (S, M & L) to fit all patients. The tong arrives sterile and ready to use in case of an emergency.

How Supplied

The STERILE set includes a carbon-graphite tong, two titanium spring-loaded pressure pins and a titanium traction hook.

Our Commitment...Customer Satisfaction

PMT® is committed to providing you the highest quality products and the finest customer service. Local sales representatives are available to serve you across the USA.

To properly size a patient: Measure a straight line 1cm above the ears, straight across the head. Refer to table below for tong sizes.

Model	Size	Description
1210-2	Small 4"-6" (10-15cm)	Supplied with 2 1/2" pins
1210-4	Medium 4 3/4"-7" (12-18cm)	Supplied with 3" pins
1210	Large 5 1/4"-7 1/2" (13-19cm)	Supplied with 2 1/2" pins


GraphComp™ 1210 Traction Tong



Features

Benefits

- | | |
|--|---|
| 1. Carbon-graphite tong with titanium pins | MRI & CT compatible |
| 2. Available in three sizes: S, M & L | Fits ALL patients |
| 3. Tong arrives sterile | Ready to use immediately |
| 4. Self-pressure pins | Do not require a torque driver or adjustment. No additional equipment needed. |

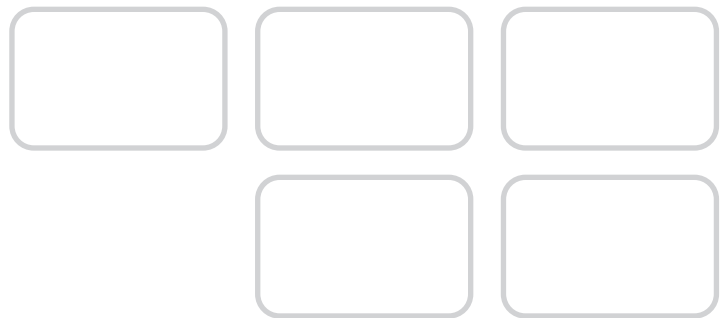

PMT® and Mankind® are registered trademarks and GraphComp™ is a trademark of PMT® Corporation. © 2008 PMT® Corporation. All rights reserved.